

CHLA-122 Cell Line Data Sheet

Cell Line Name: CHLA-122

Disease: Neuroblastoma
Phase of Therapy: Diagnosis
Prior Therapy: None
Disease Stage: 4
Source of Culture: Bone marrow
Primary Tumor Site:
Date Established: November 1992

MYCN Status:
p53 Status: Functional
TH expression: Expressed
Karyotype:
Modal No:

LC90 (DIMSCAN*):
 *see reference 4

<u>CBDCA (µg/ml)</u>	<u>ETOP (µg/ml)</u>	<u>L-PAM (µg/ml)</u>
0.1	0.1	1.6

CBDCA, carboplatin; ETOP, etoposide; L-PAM, melphalan

Gender: Female
Age: 24 months
Race:

STR Profile: On file at <http://strdb.cogcell.com>

Media Formulation: Cells are grown in a base medium of Iscove's Modified Dulbecco's Medium plus the following supplements (to a final concentration): 20% Fetal Bovine Serum, 4mM L-Glutamine, 1X ITS (5 µg/mL insulin, 5 µg/mL transferrin, 5 ng/mL selenous acid)

Growth Conditions: 37.0°C, 5% CO₂, 20% O₂
Subculturing: Next page or Protocols section at <http://www.cogcell.org/protocols.php>
Doubling Time: 72 hours
Morphology: Small, round cells
Growth Properties: Suspended, grow mostly in tight clumps

Notes: COGcell.org has a post-treatment cell line available from this same patient (CHLA-136)

All COG Repository cell lines are antibiotic-free, mycoplasma-free, and cryopreserved in 50% FBS / 7.5% DMSO. Each vial label contains the cell line name, passage number, total viable cell count (usually 5-10e6), the overall cell viability, and date frozen. All cell lines are validated with original patient sample by STR analysis.

CHLA-122 Cell Line Data Sheet

Cell Line Name: CHLA-122

GROWING CELLS FROM FROZEN CULTURES

To ensure that cells are recovered at the highest viability possible, immediately process cells and place them in culture. If the cells cannot be processed immediately, store the frozen vial in liquid nitrogen vapor (-196°C) upon receipt. Storing cells at -80°C will decrease the viability of the cells.

1. Thaw the vial in a water bath (37°C) for approximately two minutes. To prevent contamination, place the vial in a small ziplock bag or be sure to keep the O-ring and cap out of the water. Thawing should take no longer than 2-3 minutes.
2. While the vial is thawing, place 9 mL of the recommended complete cell culture medium in a centrifuge tube.
3. After the vial has thawed, record the label contents prior to processing. Spray the outside of the vial with 70% ethanol. Remove the contents of the vial and place it in the centrifuge tube containing the complete medium. Spin the tube at 500 x g for 5-7 minutes.
4. After spinning, decant the supernatant and resuspend the pellet in 13 mL of the complete medium. Transfer the cell suspension to a 75 cm² culture flask. Incubate the cells at 37°C in 5% CO₂ and room oxygen tension (~20-21% O₂).

MEDIUM RENEWAL

Protocols for medium renewal and subculturing are for a 75 cm² flask. If using different sized flasks or other cell culture vessels, adjust volumes accordingly. Cells must be renewed with new medium or subcultured at least every 3-5 days. If cells are not ≥80% confluent after 3-5 days, replace the medium with fresh medium and incubate until they are ready to be subcultured.

1. Remove supernatant from flask and place it in a 15 mL centrifuge tube. Spin the medium at 500 x g for 5 to 7 minutes.
2. After spinning, resuspend the pellet in 13 mL of complete medium and place the suspension back into the 75 cm² flask. Be sure to slowly place the medium back into the flask as to not disturb attached or loosely attached cells.

SUBCULTURING CELLS

1. Add 10-13 mL of Puck's EDTA (see formulation below) to the flask containing the cells. (CHLA-122 cells grow in tight, suspended clumps and Puck's facilitates to break up the floating spheres). Place the flask in the incubator for 10-15 minutes.
2. After incubation, gently pipette the cells up and down multiple times to help break up the clumps. (If a cell count is required, cells can be counted after this step).
3. Place the cell suspension in a 50 mL centrifuge tube and spin at 500 x g for 5-7 minutes.
4. Decant supernatant and resuspend cells in fresh complete medium and add the appropriate aliquots of the cell suspension to new culture vessels. Seed cells at a ratio of 1:3 to 1:6.
5. Incubate the cells at 37°C in 5% CO₂ and 20% O₂.

PRESERVING CELLS

Cells should be frozen in a mixture of 50% FBS, 7.5% DMSO, and base medium (with no supplements) and stored in liquid nitrogen vapor.

PUCK'S EDTA

Puck's EDTA can be made using the following cell culture grade chemicals at defined concentrations: 140 mM NaCl, 5 mM KCl, 5.5 mM Glucose, 4 mM NaHCO₃, 13 μM Phenol Red, 0.8 mM EDTA, and 9 mM HEPES.

The solution should be sterile filtered (0.22 μm) before use. The optimal pH for the solution is a range between 7.2-7.3.

CHLA-122 Cell Line Data Sheet

Cell Line Name: CHLA-122


References:

1. Keshelava N, Zuo JJ, Chen P, Waidyaratne SN, Luna MC, Gomer CJ, Triche TJ, Reynolds CP: Loss of p53 function confers high-level multi-drug resistance in neuroblastoma cell lines. *Cancer Res.* 61:6185-6193, 2001. PubMed ID: [11507071](#)
2. Keshelava N, Davicioni E, Wan Z, Ji L, Sposto R, Triche TJ, Reynolds CP. Histone Deacetylase 1 Gene Expression and Sensitization of Multidrug-Resistant Neuroblastoma Cell Lines to Cytotoxic Agents by Depsipeptide. *J Natl Cancer I.* 99: 1107-19, 2007. PubMed ID: [17623797](#)
3. Harned TM, Kalous O, Neuwelt A, Loera J, Ji L, Iovine P, Sposto R, Neuwelt EA, Reynolds CP: Sodium Thiosulfate (STS) administered six hours after cisplatin does not compromise anti-neuroblastoma activity. *Clin Cancer Res.* 14:533-540, 2008. PubMed ID: [18223229](#)
4. Kang MH, Smith MA, Morton CL, Keshlava N, Houghton PJ, Reynolds CP. National Cancer Institute Pediatric Preclinical Testing Program: Model Description for In Vitro Cytotoxicity Testing. *Pediatr Blood Cancer.* 56: 239-249, 2011. PubMed ID: [20922763](#) (www.PPTPinvitro.org)


CHLA-122 Cell Line Data Sheet

Cell Line Name: CHLA-122


Low confluency (10x magnification)

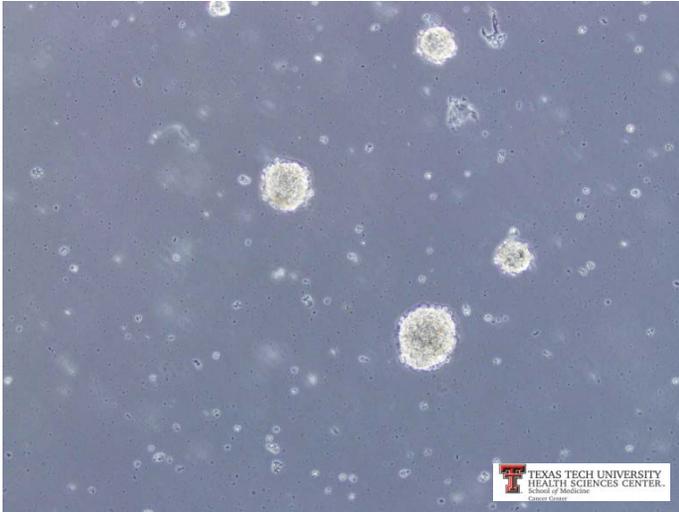

High confluency (10x magnification)

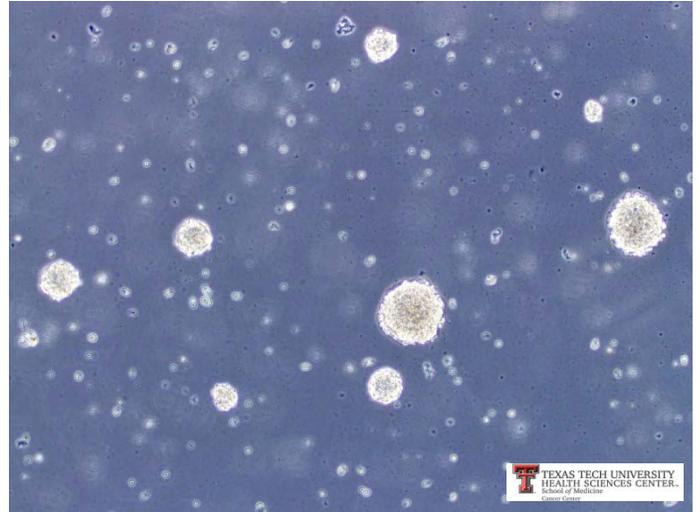

Low confluency (20x magnification)

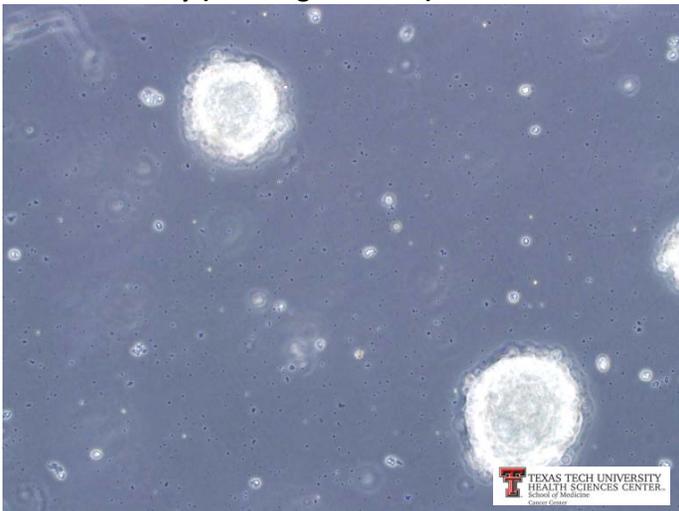

High confluency (20x magnification)

